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Version 01/11

Hochpräzise Taschenwaage PS-20

Best.-Nr. 12 41 35

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High-precision pocket scale PS-20

Item-No. 12 41 35

Use in accordance with manufacturer's instructions

The pocket scale is intended to be used for accurately weighing objects up to a total weight of 20 grams. A tare, additive weighing and adjustment function are built in.

The scale is only designed for short, direct weighing. Dynamic weighing (e.g. pouring material) is not feasible and will result in erroneous measurements.

The pocket scale is powered by batteries. The voltage supply is delivered by 4 micro batteries (AAA). These are not enclosed.

Do not operate the scale where there is a risk of explosions. Do not use the scale for humans or animals, for medical purposes, or in commercial trade for determining prices.

Any other use than described above is not permitted, it could cause damage to the product and personal injury. No part of the product may be modified or converted!

The safety instructions are to be observed without fail!

Safety Instructions

Please read all of the operating instructions before using the product for the first time; these instructions contain important information about the correct operation.

The warranty will be void in the event of damage caused by failure to observe these safety instructions! We do not assume any liability for any resulting damage!

Nor do we assume any liability for material and personal damage caused by improper use or non-compliance with the safety instructions! The warranty will be void in such cases.

This device left the manufacture's factory in a safe and perfect condition.

We kindly request the user to observe the safety instructions and warnings contained in the enclosed operating instructions so this condition is maintained and to ensure safe operation. Please pay attention to the following symbols:

An exclamation mark in a triangle points out important information in these operating instructions that have to be taken notice of.

This device conforms to CE and meets the necessary national and European regulative criteria.

The „hand“ symbol is used to indicate where specific hints and information on handling are given.

The unauthorized conversion and/or modification of the product is inadmissible for reasons of safety and approval (CE).

Avoid strong mechanical stresses, such as impacts, vibrations or drops.

Avoid strong electromagnetic fields (e.g. large motors, etc.), draughts or thermal fluctuations, as well as humidity. These may lead to measurement variations and erroneous measurements.

Keep the unit out of the reach of children. It is not a toy.

Do not leave packaging materials unattended. They could become a dangerous playthings for children.

On industrial sites, the accident prevention regulations of the association of the industrial workers' society for electrical equipment and utilities must be followed.

Never switch the device on immediately after it has been taken from a cold into a warm environment. Condensation that forms might destroy your device. Allow the device to reach room temperature before switching it on.


You should also heed the safety instructions in each chapter of these instructions.

Package Contents

Pocket scale
Reference weight 20 grams
Case
Operating Instructions

Description of the parts

- 1 Flap lid
- 2 Weighing plate
- 3 Display
 - A Battery replacement display
 - B Sign (minus)
 - C Active tare function
 - D Zero indication, measuring is possible
 - E Weight indication
 - F Indication of the weight unit
- 4 MODE button to select the weight unit
- 5 ADD UP button for adding up
- 6 TARE-button for taring and adjusting
- 7 ON-OFF button, to switch on and off
- 8 Battery compartment on the rear

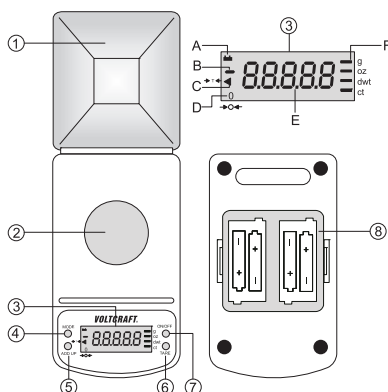Getting started

Inserting/Replacing the Batteries

The pocket scale is powered by batteries. If you use the product for the first time, or if the battery replacement symbol (A) appears on the display, you have to insert new batteries in the battery compartment.

To insert and replace the batteries, proceed as follows:

Close the flap lid (1).

Remove the lid of the battery case on the rear side (8).

Remove the used batteries from the battery compartment.

Insert four new micro batteries (AAA) with the polarity aligned. Take note of the polarity information in the battery compartment.

Close the battery compartment carefully.

Do not mix old and new batteries. Always replace all the batteries at the same time, otherwise the batteries discharge each other. Always use batteries of the same type and brand.

You can acquire suitable micro batteries under following order number:
order no. 65 23 47 (package of 4, please order 1x).

Do not leave flat batteries in the device. Even batteries protected against leakage can corrode and thus release chemicals that could be detrimental to your health or destroy the device.

Remove the batteries, if the appliance is not to be used for long while, to prevent leaking.

Leaking or damaged batteries might cause severe burns if they come into contact with the skin. Use appropriate protective gloves in such case.

Do not leave batteries lying around carelessly. They might be swallowed by children or pets. If swallowed, consult a doctor immediately.

Make sure that the batteries are not short-circuited. Do not throw batteries into a fire.

Batteries (non-rechargeable) must not be recharged. Danger of explosion.

Weighing

The maximum weighing range of 20 g must not be exceeded. Avoid thermal fluctuations during weighing; this may lead to erroneous measurements.

Preparation

Place the pocket scale on a stable, level surface and open the flap lid.

Press the „ON/OFF“ button (7), to switch on the scale. After a brief testing phase (all digits are displayed), the scale performs the zero setting (display „0.00“ or „0.000“). The scale is now ready for use. Avoid vibrations and loading of the weighing plate during this phase. When the check has been completed successfully, a circle appears in the lower left corner (-> o <-) indicating stable measuring conditions, as well as the reading „0.00“ or „0.000“.

If the automatic zero setting has failed, remove any weights that might be placed on the weighing plate and make sure it is on an even surface. Avoid vibrations and draught. Turn the device off and then on again.

If these measures do not bring the desired result, the zero setting can be carried out manually using the „TARE“ button (6).

Press the „ON/OFF“ button to turn off.

Selecting the Weight Unit

The pocket scale provides the option of weighing with 4 different weight units.

Gram (g), ounce (oz), pennyweight (dwt) or carat (ct). Press the „MODE“ button (4), when the scale is switched on. Each press switches the weight unit. The selected unit is shown by a line, on the right side of the display. The last selected unit remains selected after switching the device off.

Weighing

Switch on the scale and, depending on the unit, wait until „0.000“ or „0.000“ is displayed.

Open the flap lid (1).

Place the object to be weighed carefully on the weighing plate (2).

Close the flap lid to avoid erroneous measurements due to slight air movements. Make sure the weighing material does not touch the scale housing or the flap lid.

Wait for approx. 5 seconds, until the reading is constant. The measurement result is displayed.

Press the „ON/OFF“ button to turn off.

When the measuring range is exceeded, „- - - -“ is shown. In this case, immediately remove the item to be weighed from the pocket scale, as it may be damaged otherwise. Remove the object from the scale, immediately after measuring, to avoid damage to the scale.

Tare function

To weigh liquids or bulk goods, weighing vessels can be used. To correct for the weight of the weighing vessel during measuring, the tare function can be used.

The total weight of weighing vessel and weighing object must not exceed the max. weight of 20 g. Thus, use light plastic or paper vessels.

Proceed as follows to weigh:

Switch on the scale, and wait until the zero indication is displayed.
Place the empty weighing vessel on the weighing plate (2) and wait until the displayed value is constant.
Then briefly press the "TARE" button (6). The displayed value is set to zero. Thus, the weight of the weighing vessel is not included in the subsequent measurement. The active tare function is indicated by a small triangle in the left display frame (-> T <-).
Remove the weighing vessel and fill with the material to be weighed.
Place the filled vessel carefully on to the weighing plate (2). Dynamic measurements, i.e. by filling with material is not feasible, as this leads to erroneous measurements.
Wait for approx. 5 seconds, until the reading is constant. The measurement result is displayed.
The Tare function can be switched off by pressing the „TARE“ button again.
Press the „ON/OFF“ button to turn off.

Additive weighing function

The additive weighing function enables the adding of subsequent measurements. With this function individual weights of max. 20 g each, up to 99.999 g can be added together. This makes it possible to create the precise mixture of different substances.
Proceed as follows to weigh:
Switch on the scale, and wait until the zero indication is displayed.
Open the flap lid (1).
Place the first object to be weighed carefully on the weighing plate (2). Wait for approx. 5 seconds, until the reading is constant. Press and hold the „ADD UP“ button (5) for approx. 2 s. The measured value will be saved. Remove the object from the scale. Wait until zero is shown on the display.
Place the next object to be weighed carefully on the weighing plate (2). Wait for approx. 5 seconds, until the reading is constant. Press and hold the „ADD UP“ button (5) for approx. 2 s. The measured value will be added up, displayed and saved. Remove the object from the scale. Wait until zero is shown on the display.
Proceed with this step for your other additive measurements.
The total weight can be read out for approx. 2 seconds after the last additive measurement. If the total weight of 99.999 is exceeded, the message „Err“ is shown on the display.
To delete the total weight saved, press and hold the „ADD UP“ until „CLF“ appears on the display, while the scale and zero indication are switched on.
Press the „ON/OFF“ button to turn off.

Adjustment

The pocket scale should be re-adjusted after a while, to avoid measuring inaccuracies. For the adjustment, a reference weight of 20 g is included in the delivery.
To carry out the adjustment, proceed as follows:
Switch the scale off and place it on a level and stable surface.
Press, and hold, the „TARE“ button for approx. 2 seconds, until a numerical value appears on the display. An exact value is not relevant.
Wait until the value has stabilised and press the „TARE“ button again.
„ZErø“ is shown on the display, and then the calibration value „20“.
Place the reference weight of 20 g on the scale plate and then press the „TARE“ button.
In the display „PASS“ will appear. Remove the reference weight from the scale and switch the scale off. The scale is readjusted.
 Treat the reference weight with great care. If possible, only touch it with gloves. Any soiling or damage leads to weight variations, and consequently to inaccuracies during the adjustment.

Automatic cut-off

To save batteries, and to increase the operating time, the pocket scale switches off automatically after approx. 3 - 4 minutes, if no keystroke is performed and no weight is measured.
Press the „ON/OFF“ button, to switch the scale on again.

Maintenance and Cleaning

Apart from occasional cleaning and battery replacement, the device requires no servicing.
Use a soft, dry, antistatic and lint-free cloth for cleaning. Do not use abrasive or chemical cleaners or petrol.

Disposal

Electronic products are raw material and do not belong in the household waste. When the device has reached the end of its service life, please dispose of it, according to the current statutory requirements, at your local collecting site. Disposing of the device in the household waste is prohibited.

Disposal of Used Batteries

As a consumer, you are required (**Battery Ordinance**) to responsibly dispose of all used batteries and rechargeable batteries; **it is forbidden to throw them away with the normal household waste!**

You can dispose of your used batteries/rechargeable batteries free of charge at your community's collection point or any place where batteries/rechargeable batteries are sold!
You thus fulfil the legal requirements and make your contribution to protecting the environment!

Technical Data

Weighing range max.	20 g / 0,705 oz / 12,86 dwt / 100 ct
Definition	0.001 g <10 g, 0.002 g >10 g
Accuracy	+/- (0.2 % + lowest indicated digit)
Overload display	-----
Weight units	g, oz, dwt, ct
Measuring sensor	Wear-free strain gauge
Power supply	4 micro batteries (AAA)
Automatic cut-off	approx. 3 - 4 minutes
Operating conditions	18 to 25 °C, max. 85% rel. air humidity, non-condensing
Dimensions (L x W x H)	113 x 65 x 18 mm
Weight (scale)	approx. 78 g

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© Copyright 2011 by Voltcraft*